
Panasonic®

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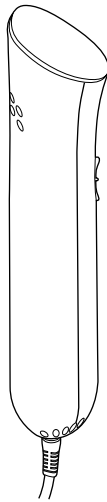Panasonic®

Operating Instructions (Household) Facial Cool Putter

Model No. EH-SQ10

Contents

Safety precautions	16
Parts Identification	22
Notices for use.....	23
How to use.....	24
Maintenance of facial cool putter	26
Troubleshooting	26
Specifications	27

Thank you for purchasing this Panasonic product.
Before operating this unit, please read these instructions completely and save them for future use.

English

Safety precautions

To reduce the risk of malfunction, electric shock, injury, fire, loss of life, and damage to equipment or property, always observe the following safety precautions.

Explanation of symbols

The following symbols are used to classify and describe the level of hazard, injury, and property damage caused when the notation is disregarded and improper use is performed.

WARNING

Denotes a potential hazard that could result in serious injury or death.

CAUTION

Denotes a hazard that could result in minor injury.

The following symbols are used to classify and describe the type of instructions to be observed.

This symbol is used to alert users to a specific operating procedure that must not be performed.

This symbol is used to alert users to a specific operating procedure that must be followed in order to operate the unit safely.

WARNING

Women in pregnancy and menstruation do not use this product.

- Doing so may cause skin problem due to hormonal imbalance.

Do not apply the unit to those persons who are insensitive to temperature, or infants.

- Doing so may cause frostbite, electric shock, or injury.

WARNING

Do not store within the reach of children or infants. Do not let them use this product.

- Putting the protective cover or any other part in the mouth may cause an accident or injury.

Do not use in a way that exceeds the rating of the household outlet or the wiring.

- Exceeding the rating by connecting too many plugs to one household outlet may cause fire due to overheating.

Do not connect or disconnect the AC adaptor to a household outlet with a wet hand.

- Doing so may cause electric shock or fire due to a short circuit.

Do not use anything other than the supplied AC adaptor.

Also, do not charge any other product with the supplied AC adaptor.

- Doing so may cause electric shock, burns, or fire due to a short circuit.

Never use the appliance if the AC adaptor is damaged or if the power plug fits loosely in a household outlet.

- Doing so may cause electric shock or fire due to a short circuit.

Do not damage, break, modify, or forcibly bend, pull, or twist the cord.

Also, do not place anything heavy on or pinch the cord.

- Doing so may cause electric shock or fire due to a short circuit.

WARNING

Do not repair, disassemble or modify the product by yourself.

- Doing so may cause fire or injury due to abnormal action. Contact a service center authorized by Panasonic for repair when necessary.

Do not use this product near other inflammable materials (such as alcohol, gasoline, diluent, spray, hair care articles, or nail polish remover, etc.).

- Doing so may cause explosion or fire.

Do not use this product in a bathroom or shower room.

- Doing so may cause electric shock or fire due to a short circuit.

Do not immerse this product or AC adaptor in water or wash with water.

Do not place or store the unit in places adjacent to a sink, water, or other highly humid places containing water.

Keep the appliance dry.

- Failure to do so may cause electric shock or fire.

The following persons should consult a qualified physician prior to using this appliance.

- Patients with acute illness.
 - Patients with a malignant tumor.
 - Persons with allergies.
 - Persons with sensitive skin.
 - Persons with skin problems, such as dermatitis severe sunburn, etc.
 - Persons with skin injuries, eczema, swelling, etc.
- Failure to do so may result in problems with the skin or body.

WARNING

Immediately stop using this product if you feel pain or discomfort.

- Failure to do so may result in problems with skin or body.

Make sure to use only 100-240 V ~ voltage.

Fully insert the AC adaptor or appliance plug.

- Failure to do so may cause electric shock, injury, or fire.

After use, always set the main switch to position "0", and unplug the AC adaptor from the household outlet.

- Failure to do so may cause burns, electric shock, or fire.

Unplug the AC adaptor from the household outlet if the product is not in use.

- Failure to do so may cause electric shock or fire due to insulation failure.

Regularly clean the power plug and the appliance plug to prevent dust from accumulating.

- Failure to do so may cause fire due to insulation failure caused by humidity.

Immediately stop using and remove the AC adaptor if there is an abnormality or failure.

- Using this appliance in such a condition may cause electric shock, injury, or fire.

<Abnormality or failure cases>

- The main unit, AC adaptor, or cord is deformed or abnormally hot.
- The main unit, AC adaptor, or cord smells of burning.
- There is abnormal sound during use or charging of the main unit, AC adaptor, or cord.
- Immediately contact a customer service center.

WARNING

The supply cord cannot be replaced. If the cord is damaged, the AC adaptor should be disposed of.

- Failure to do so may cause an accident or injury.

Properly dispose of the packaging of the appliance, keep the packaging away from infants.

- Failure to do so may cause suffocation of infants.

This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

Children should be supervised to ensure that they do not play with the appliance.

- Failure to do so may cause an accident or injury.

CAUTION

Do not allow pins or other small objects to adhere to the AC adaptor, the appliance plug, or the appliance socket.

- Doing so may cause electric shock or fire due to a short circuit.

Do not drop or subject to shock.

- Doing so may cause injury.

Do not wrap the cord around the AC adaptor when storing.

- Doing so may cause cord wire breakage, electric shock, or fire due to a short circuit.

Do not use this product below 15 °C.

- Doing so may cause skin frostbite should your skin come in contact with the cold plate.

CAUTION

Check if the cold plate is damaged or deformed before use.

- Failure to do so may cause injury.

Only use this appliance for face are.

- Failure to do so may cause bodily injury or discomfort.

Hold the AC adaptor or the appliance plug instead of only the cord while unplugging.


- Doing so may cause cord wire breakage, burns, electric shock, or fire due to a short circuit.


This product should be cleaned after use, to prevent grease or other residues accumulating.

- Failure to do so may cause injury to your skin or body.

Parts Identification

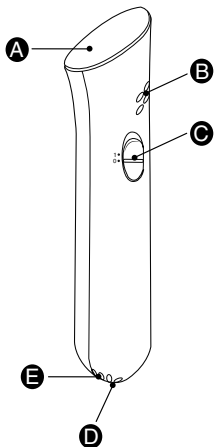
- A** Cold plate
- B** Air inlet
- C** Main switch [0/1]
1:ON
0:OFF
- D** Appliance socket
- E** Air outlet

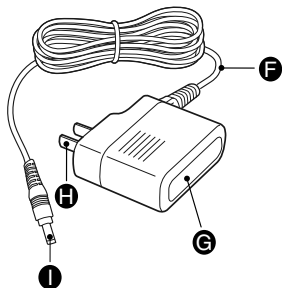- F** Cord
- G** AC adaptor
- H** Power plug
- I** Appliance plug

Notices for use

- At the beginning, it is advised to use this product no more than 2 - 3 times in a week.
You can use it once a day (about 2 minutes) after a month.
Adjust the frequency and duration depending on your skin condition.
- Skin flush or pruritus may occur during use.
Immediately see a doctor if the symptom is not mitigated after you stop using this product.
- Do not use this product when wearing makeup.
- Use this product only in the range of 15 ~ 40 °C.
Cold injury may occur below 15 °C. Product failure is more likely to occur at higher temperature.
- Do not apply this product to the same part of your body for a long time.
- Do not block the air inlet or the air outlet (the shaded portion in the illustration).
Doing so may cause not cold on the cold plate.

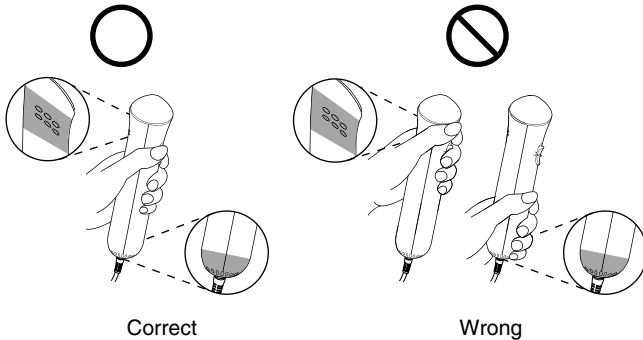How to use

- Check if the surface of the cold plate is damaged or dirty before use.
- 1. Plug the appliance plug into the appliance socket and insert the AC adaptor into a household outlet.**
 - 2. Set the main switch to “1”.**
Operation sound should be audible when power is turned on.
The cold plate will become cold after about 30 seconds.
 - 3. Gently touch your face with the cold plate.**
It will take about 1 minute to finish the cold compress on your whole face.

■ Tips

- Move over your face with a gentle patting motion.
- Move from inside to outside of your face.
- Actions on your left face and right face should be symmetrical.

Application diagram


Use the entire cold plate

For wide parts of your face, such as the cheeks and forehead

Use the tip of the cold plate

For narrow parts of your face, such as side of the nose

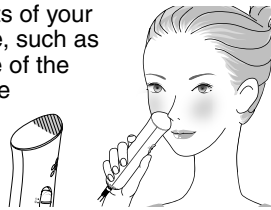
- 4. Set the main switch to “0” and unplug the AC adaptor from the household outlet after cold compressing.**
The cold plate may be temporarily warm, but will restore to normal temperature after power-off.

Recommended method to use with facial steamer


■ Notice for use

- Do not directly blow this product by steam.

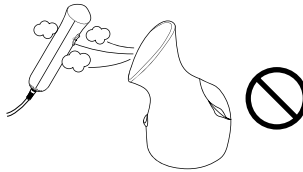■ Recommended use mode (take the facial steamer running for 6 minutes as an example)

● As a final step of moisturizing

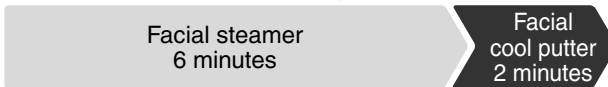(Set the main switch to "1" one minute before using the facial cool putter)

● On pore care

■ After use

- Use moisturizing skin care products to take care of your face after using this appliance. If you do not take care of your skin, water will evaporate from it, making it dry.

Maintenance of facial cool putter

Do not use paint thinner, benzene, alcohol, etc., to clean the unit. Doing so may cause failure, cracking or discoloration of the main body. Wipe with a cloth soaked in soapy water and wrung tightly.

Troubleshooting

Problem	Possible cause	Action
Not cold	▶ Main switch is off	▶ Set the main switch to "1"
	▶ The air inlet or air outlet is blocked	▶ Do not block the air inlet or air outlet while holding the appliance
	▶ The cord shows signs of breakage ▶ The cord is disconnected	▶ Contact a customer service center

If the problems still cannot be solved, contact a customer service center.

Specifications

Product name		Facial cool putter
Model No.		EH-SQ10
Rated input of AC adaptor		100-240 V ~ 50-60 Hz 180 mA
Rated output of AC adaptor		5.4 V \equiv 1.2 A
Place of origin		China
Date of production		See engraving on the back of the commodity.
Execution standard	Chinese standard	GB4706.1 GB4706.15
	Company standard	Q/SXWBMJ 15

This product is intended for household use only.

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Panasonic Wanbao Appliances Beauty and Living (Guangzhou) Co., Ltd.