Panasonic

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English

Panasonic

Operating Instructions
(Household) IPL Hair Removal System
Model No. ES-WH93

Contents

Safety precautions 24	For healthy skin30
Intended use29	<u>-</u>
Parts identification30	Cleaning3
Treatment areas31	Q&A38
Guide for treatment frequency, time31	Troubleshooting40
Charging32	Battery life40
Before use33	Removing the built-in rechargeable
When using the appliance for the	batteries4
first time34	Specifications42
How to use34	

Thank you for purchasing this Panasonic product.

Before operating this unit, please read these instructions completely and save them for future use.

Safety precautions

To reduce the risk of injury, loss of life, electric shock, fire, and damage to property, always observe the following safety precautions.

Explanation of symbols

The following symbols are used to classify and describe the level of hazard, injury, and property damage caused when the denotation is disregarded and improper use is performed.

⚠ DANGER

Denotes a potential hazard that will result in serious injury or death.

⚠ WARNING

Denotes a potential hazard that could result in serious injury or death.

⚠ CAUTION

Denotes a hazard that could result in minor injury or property damage.

The following symbols are used to classify and describe the type of instructions to be observed.


This symbol is used to alert users to a specific operating procedure that must not be performed.

This symbol is used to alert users to a specific operating procedure that must be followed in order to operate the unit safely.

⚠ WARNING

► To prevent skin or health problems

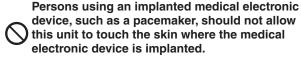- Doing so may cause an accident or illness.

Do not use during pregnancy or menstruation.

Doing so may cause skin problems due to unstable hormone balance.

Do not use on tattoos.

- Doing so may cause keloids by burns or redness or other rashes on skin.

Do not shine directly into the eyes.

- Doing so may cause injury to the eyes.

Only use on underarms, arms, legs, bikini-line or face (when using the face attachment).

- Failure to do so may cause burns or redness or other rashes on skin

► Power supply

Do not connect or disconnect the adaptor to a household outlet with a wet hand.

- Doing so may cause electric shock or injury.

WARNING

Do not damage or modify, or forcefully bend, pull, or twist the cord. Also, do not place anything heavy on or pinch the cord.

Never use the appliance if the AC adaptor is damaged or if the power plug fits loosely in a household outlet.

 Doing so may cause electric shock or fire due to a short circuit.

Do not use anything other than the supplied AC adaptor for any purpose.

Also, do not use any other product with the supplied AC adaptor. (See page 30.)


 Doing so may cause burn or fire due to a short circuit.

Do not use in a way exceeding the rating of the household outlet or the wiring.

 Exceeding the rating by connecting too many plugs to one household outlet may cause fire due to overheating.

Always ensure the appliance is operated on an electric power source matched to the rated voltage indicated on the AC adaptor.
Fully insert the adaptor or appliance plug.

- Failure to do so may cause fire or electric shock.

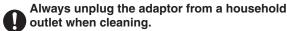- Failure to do so may cause electric shock or injury.

Regularly clean the power plug and the appliance plug to prevent dust from accumulating.

- Failure to do so may cause fire due to insulation failure caused by humidity.

Disconnect the adaptor and wipe with a dry cloth.

▶This product

This product has a built-in rechargeable battery.

Do not throw into flame or heat. Do not charge,
use, or leave in high temperatures.

Doing so may cause overheating, ignition, or explosion.

Do not store within the reach of children. Do not let them use it.

- Doing so may cause an accident or injury.

Do not modify or repair.

Ontact an authorized service center for repair (battery change etc).

Never disassemble except when disposing of the product.

- Doing so may cause fire, electric shock, or injury.

On not use if the air inlet is blocked.

- Doing so may cause fire or burns.

WARNING

► Main body and AC adaptor

Do not store the main body and AC adaptor in a bathroom or location with high humidity. Also, do not place in a location where water may spill (such as top of a vanity).

- Doing so may cause electric shock or fire.

Do not immerse the main body and AC adaptor in water or wash it with water.

Do not use the main body and AC adaptor in the bathroom or location with high humidity.

 Doing so may cause electric shock or fire due to a short circuit.

Keep the appliance dry.

► In case of an abnormality or malfunction

Immediately stop using and remove the adaptor if there is an abnormality or malfunction.

- Using it in such conditions may cause fire, electric shock, or injury.
- <Abnormality or malfunction cases>
- The main unit, adaptor or cord is deformed or abnormally hot.
- •The main unit, adaptor or cord smells of burning.
- There is abnormal sound during use or charging of the main unit, adaptor or cord.
- Immediately request inspection or repair at an authorized service center.

► Take care of the following precautions as well

This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction

been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

Children should be supervised to ensure that they do not play with the appliance.

- Failure to do so may cause an accident or injury.

The supply cord cannot be replaced. If the cord is damaged, the AC adaptor should be scrapped.

- Failure to do so may cause an accident or injury.

Properly dispose the package of the appliance, keep the package away from infants.

-Failure to do so may cause suffocation of infants from happening.

ACAUTION

► To protect the skin

The following persons should not use this appliance:

- · Persons with dermatitis
- · Persons suffering from atopic dermatitis or rashes and other skin diseases
- Persons with allergic diathesis, or sensitive skin which is easily irritated by cosmetics, clothing, metals, etc.
- Persons with skin that is sensitive to sunlight

Do not use in the following body area:

- · On suntanned skin
- Dark brown or black skin areas such as moles See the sticker of "Suitable skin tones" on the back of the main body.

- Scabs, injuries, rashes, bruises, warts, pimples, skin cancer, ringworm (dermatophytosis), folliculitis, irritated skin, insect bites, etc.
- · Angioma
- Areas with cosmetic surgery

Do not use in the following cases:

· When you are not well or when you feel an abnormality in your skin or body

Do not use on the same area repeatedly. (Try using the appliance up to once every 3 days even if removal of unwanted hair is not immediately obvious.)

 Doing so may cause injury to the skin, inflammation, bleeding, increased pain or other problems or aggravated symptoms.

Do not remove the attachment except during cleaning.

- Doing so may cause injury to the skin or fingers.

Only for use by persons over the age of 18.

- Failure to do so may cause injury or inflammation of the skin, bleeding or increased pain.

Before use, check that there is no distortion, cracks or damage to the light guide and flash window.

- Failure to do so may cause injury to the skin.

•

Persons on medication should consult a dermatologist.

- Skin may become sensitive to light stimulation depending on the type of medication being taken, and may cause skin problems.

If inflammation (such as skin turning red) occurs, stop using the appliance, and consult a dermatologist if the inflammation continues even after 2 days.

- Failure to do so may cause symptoms to worsen.

Use after removing metal objects such as accessories.

- Failure to do so may cause an accident or injury.

► Power supply

Do not allow metal objects or trash to adhere to the power plug or the appliance plug.

 Doing so may cause electric shock or fire due to a short circuit.

ACAUTION

Do not wrap the cord around the adaptor when storing.

- Doing so may cause the wire in the cord to break with the load, and may cause fire due to a short circuit.

Disconnect the adaptor or the appliance plug by holding onto the adaptor or the appliance plug instead of the cord.

- Disconnecting by holding onto the cord may cause electric shock or injury.

Disconnect the adaptor from the household outlet when not charging.

- Failure to do so may cause electric shock or fire due to electrical leakage resulting from insulation deterioration.

► This product

Do not use on wet skin.

- Doing so may cause an electric shock.

Do not drop or subject to shock.

- Doing so may cause injury.

Do not remove the attachment or touch the flash window immediately after using IPL hair removal.

Doing so may cause burns as the flash window becomes hot.

► Note the following precautions

Do not use on clothes or head hair.

-Doing so may cause burns or discoloration.

Handling of the removed battery when disposing

A DANGER

The rechargeable battery is exclusively for use with this appliance. Do not use the battery with other products.

Do not charge the battery after it has been removed from the product.

- Do not throw into fire or apply heat.
- · Do not hit, disassemble, modify, or puncture with a nail.

- Do not let the positive and negative terminals of the battery get in contact with each other through metallic objects.
- Do not carry or store the battery together with metallic jewellery such as necklaces and hairpins.
- •Do not use or leave the battery where it will be exposed to high temperatures, such as under direct sunlight or near other sources of heat.
- Never peel off the tube.
- Doing so may cause overheating, ignition, or explosion.

WARNING

After removing the rechargeable battery, do not keep it within reach of children and infants.

0

-The battery will harm the body if accidentally swallowed.

If this happens, consult a physician immediately.

If the battery fluid leaks out, do not touch the battery with your bare hands.

-The battery fluid may cause blindness if it comes in contact with your eyes.

Do not rub your eyes. Wash immediately with clean water and consult a physician.

-The battery fluid may cause inflammation or injury if it comes in contact with the skin or clothes. Wash it off thoroughly with clean water and consult a physician.

Intended use

- Please observe the following before performing IPL hair removal.
 - This appliance should only be used on bare skin with no cosmetics, lotions etc.

When using the lotion or shaving agents (gels, etc.), rinse off the agent thoroughly with water and let the skin dry.


- · Be sure to shave the area (body or face). Before use, wipe away any water or sweat from the treatment area, and any hair sticking to the skin first. Hairs or shaving agents stuck to the skin can foul the light guide. It may become difficult to clean the light quide.
- Do not use alcohol, nail polish remover, or detergent (ex. liquid hand soap), etc. to clean the appliance.
 Doing so may cause failure, cracking or discoloration of the main body. Wipe with a dry cloth.
- Do not leave the appliance in an area with high temperature.

Doing so may cause malfunction.

 The following symbol indicates that a specific detachable power supply unit is required for connecting the electrical appliance to the supply mains. The type reference of power supply unit is marked near the symbol.

Parts identification

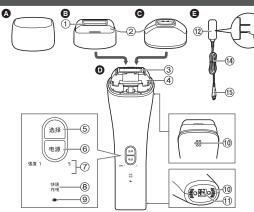- **A** Cap
- Body attachment
- **©** Face attachment
 - 1 Light guide
 - ②Flash switch
- Main body
 - ③Flash window
 - 4 Flash lamp
 - (5) Level selector switch (flashes)
 - 6 Power switch
 - ①Level indication lamp (blue) (flashes) (See pages 34 and 36.)

- ® Quick charging indication lamp (red) (See pages 32 and 37.)
- 10 Air inlet
- 11 Socket
- **⊕** AC adaptor (RE7-87)
 - 12 Adaptor
 - 13 Power plug
 - (14) Cord
 - (15) Appliance plug

Removing and mounting the body/face attachment

▶ Removing the attachment

Holding the attachment by side, remove from the main body.

► Mounting the attachment

Push the attachment until it clicks.

Treatment areas

Face

Do not use on the area, such as around the eye, eyebrow and lip, etc.

- Tie up hair or fix it with a hair band or similar. Also, please also fix short hair with a clip or similar.
- Please take care in particular of sunburn before and after using on your face.
- Do not use for the man's beard.

Guide for treatment frequency, time

*1 For bikini-line treatment, wear a swim wear, under wear etc. first to check the treatment area.

*2 Do not use in the abdomen and genitals.

Frequency

Legs/Arms

Use once every 2 weeks for the first 2 months.

Use once every 4 weeks from the 3rd month onward.

Underarms/Bikini-line/Face

Use once every week for the first 2 months.

Use once every 2 weeks from the 3rd month onward.

Underarms

Arms

Bikini-line *1,2

Legs


Guide for treatment time

Legs	Approximately 10 minutes per session for both legs (approximately 190 flashes)
Arms	Approximately 5 minutes per session for both arms (approximately 90 flashes)
Underarms	Approximately 2 minutes per session for both underarms (approximately 30 flashes)
Bikini-line	Approximately 2 minutes per session (approximately 30 flashes)
Face	Approximately 3 minutes per session (approximately 40 flashes)

- Note that, for facial skin care, it is recommended to use the appliance once every 3 days.
- Do not use on the same area repeatedly, even if removal of unwanted hair is not immediately obvious. Try using the appliance up to once every 3 days. (Treatment effects may differ depending on the area or individual body types.)
- If the appliance has not been used for a while, resume use at the same frequency as the first 2 months.
- Treatment time may differ depending on the area or amount of hair, treatment area and how the appliance is used.
- Please shave any unwanted hair first from the treatment area. Please use a separate razor or electric shaver.
- The appliance differs from permanent hair removal treatment. Unwanted hair will gradually grow back if treatment is stopped.


Charging

- Insert the appliance plug firmly into the socket at the bottom of the main body.
- **2** Firmly insert the adaptor into a household outlet.

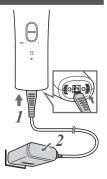- 3 Check that the charging indication/notification lamp glows.
 - If the remaining capacity of the battery is low, the appliance will start from quick charging. (Quick charging is completed after approx. 10 minutes.) Depending on the remaining capacity of the battery, the appliance may start from normal charging rather than quick charging.

The quick charging time may differ with the charging capacity.

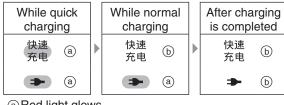- a Red light glows.
- b Red light goes off.
- Charging is completed in about 3 hours.
- The charging time may differ with the charging capacity.
- <To check if charging is completed>
 Press the power switch with the AC adaptor connected. The charging indication/notification lamp glows, and goes off after approximately 5 seconds.
- 4

Remove the adaptor from the household outlet, and remove the appliance plug from the main body. (for safety and to reduce energy usage)

•There will be no problems with the main body even if it is charged for 3 hours or more.

Notes

- The charging time may differ when the appliance is being used for the first time, or if it has not been used for half a year or more. The lamp on the main body may not glow for several minutes after charging starts, but will glow if charging continues.
- The appliance can be used for approximately 600 flashes with a full single charge. (The appliance can be used for approximately 40 flashes with quick charging.)

- The number of flashes may differ depending on the ambient temperature and how the appliance is used.
- Recommended ambient temperature for charging is 15 – 35 °C. Charging may not be possible outside of the recommended temperature.
- If there is crackling noise on radios or other devices while charging, charge the appliance using a different household outlet.
- Charging the appliance every time it is used will not affect the life of the rechargeable battery.

Before use

- Do not use the appliance in dark places.
 Please close your eyes when flashing your face with the face attachment.
- Doing so may cause temporary irritation in the eyes.
- •This appliance should only be used on bare skin with no cosmetics, lotions etc.
- When using the lotion or shaving agents (gels, etc.), rinse off the agent thoroughly with water and let the skin dry.
- Be sure to shave the area (body or face).
 Before use, wipe away any water or sweat from the treatment area, and any hair sticking to the skin first.
 Hairs or shaving agents stuck to the skin can foul the light guide. It may become difficult to clean the light guide.

When using the appliance for the first time

Check that the output level suits your skin.

Treated skin may temporarily turn red or feel painful depending on the body type or condition of the person using the appliance.

Set the output level to 1 (强度1) and flash the treatment area once.

After 2 days

If there are no skin abnormalities

Set the output level to 1 (强度1) and start treatment on the area that you want.

• Guide for treatment frequency, time (See page 31.)

If there are no skin abnormalities 1 week

To increase the output level for treatment, flash once on the treatment area, and check your skin condition for around 2 days before further treatment. Increase the output level one step at a time.

• If there is intense pain or redness after 2 or more days, stop using the appliance and consult a dermatologist.

<Output levels> Greater treatment effects 强度5 ○ ○ ○ ○ 强度3 ● ● ● Minimal pain or irritation

How to use

- The power will not turn on even when pressing the switch with the AC adaptor connected to the main body.
- The appliance may not operate in an ambient temperature lower than approximately 5 °C.
- Remove the cap.

- Attach the attachment to suit the area where you want to use it.
- Press and hold (for approximately 3 seconds) to turn the power on.

• The power will turn off automatically if the appliance is not used for approximately 1 minute.


- Press (选择) to select the output level.
 - •Level 1 is selected when the power is turned on.
 - The level cycles through $1 \rightarrow 2 \rightarrow 3 \rightarrow 4 \rightarrow 5 \rightarrow 1...$ each time the level selector switch is pressed.
 - •The level indication lamp blinks, and then turns on. (A beep sounds, indicating that preparation for flash is complete.)

When the level indication lamp glows, apply the flash switch closely against the skin and keep it pressed until flashes. When the appliance beeps, press the head perpendicularly against

the skin.


- Please do not press too strongly.
- Use the appliance after checking that there are no other people nearby.
- · Light might leak from the air inlet on the back of the main body when flashing the treatment area; however, there will be no problem.
- Points for IPL hair removal (See this page.)
- · When the appliance beeps, preparation for flash is complete.

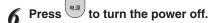Points for proper IPL hair removal

Point 1

Apply the flash switch perpendicularly against the skin and keep it pressed against the skin.

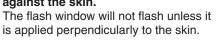
Point 2

The size of the area is a quide for the IPL hair removal range of each flash.

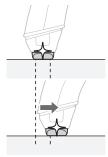- (a) Body attachment
- (b) Face attachment
- · When the appliance beeps, preparation for flash is complete.

Continuous flashes

After one flash, slide the appliance to the next flash area while it is still pressed against the skin. Up to 10 continuous flashes at an interval of approximately 3 seconds is possible.

- •The flash window will no longer flash after 10 flashes. Remove the appliance from the skin and try again.
- Do not use the appliance continuously on the same area.
- Depending on usage conditions, the flash interval may vary.

For healthy skin

The following care is recommended after IPL hair removal.

 Cool the skin using a towel soaked in ice water.

Moisturize the skin with the lotion etc.

Lamp indication

► Normal indication

強度 1、1、1、1、1、1、1、1、1、1、1、1、1、1、1、1、1、1、1、	The level indication lamp (blue) blinks (once/second). The image on the left is when level 5 is selected.
*	Preparation for flash in progress.
强度 ①	The level indication lamp (blue) glows. The image on the left is when level 5 is selected.
→	Flashing is possible.

► Charging notification

强度 1 5	(red) blinks (once/second).
快速 充	The battery level is low. Charging is recommended.

► Malfunction display

[Abnormal charging]

The charging indication/notification lamp (red) blinks quickly (twice/second).

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Recommended ambient temperature for charging is 15 – 35 °C. Check the room temperature of the place where the appliance is going to be charged, and charge it again.

[Other malfunctions]

The level indication lamp (blue), quick charging indication lamp (red) and charging indication/notification lamp (red) blink simultaneously (once/second), the appliance beeps twice, and the power turns off after 5 seconds.

The appliance may not operate in an ambient temperature lower than approximately 5 °C. Check the room temperature of the place that the appliance is going to be used, and turn it on again.

(The appliance may not adapt to the room temperature immediately after changing the location.)

If the symptom occurs repeatedly, contact an authorized service center.

Cleaning

Always turn the main body power off and wait for 2 to 3 minutes before cleaning.

► Flash window/Light guide

1. Remove the attachment.

2. Wipe any dirt clean carefully with a dry cloth.

Flash window

Light guide

3. Attach the attachment. Insert until it clicks.

► Air inlet

Wipe any dirt clean with a dry cloth.

Is it alright to use the

appliance everyday?

Q&A Question **Answer** What type of light is Instantaneous flash from a W emitted from the flash xenon lamp. sr window? This is because of unstable Why can the appliance hormone balance when be used only by those growing up and the skin can over the age of 18? become sensitive easily. Only use the appliance on br underarms, arms, legs, ap On which areas can bikini-line and face (when pr the appliance be using face attachment). used? (Do not use the body attachment on the face or other areas.) Using the appliance on the same area everyday is not

recommended.

(See page 31.)

Please follow the guide for treatment frequency and time.

Question		Answer
/hy is there a burning mell after a flash?	•	This is because the flash applies heat to unwanted hair, and is not a sign of a malfunction. Shave any hair if it has grown too long. Also remove all trimmings remaining on the skin.
he light feels very right while using the ppliance but is this a roblem?	•	Apply the flash switch closely against the skin. (See page 35.) If the light is still a concern, wear sunglasses or other eye protection.
an the flash lamp be placed?	•	The flash lamp cannot be replaced. If the flash lamp is damaged or the like, please contact the store where you purchased the unit or a service center authorized by Panasonic for repair.

Anewor

Question

Question

The light feels hot while using the appliance but is this a problem?

The treated skin turns red while using the appliance but is this a problem?

The pain and redness do not subside after use but is this a problem?

During charging, quick charging does not start? (Quick charging indication lamp does not glow)

Answer

The sensation is different for each person but it may feel hotter if hair has not been shaved properly.

- Treated skin may temporarily turn red or feel painful depending on the body type or condition of the person using the appliance. Stop using the appliance immediately, and keep the skin cool and moist. If the symptoms subside within 2 days, use the
- appliance at a level that the symptoms do not appear. If the symptoms persist after 2 days, consult a dermatologist. (See page 34.)

Depending on the remaining capacity of the battery, normal charging starts rather than Quick charging. (Charging indication/

(Charging indication/ notification lamp only glows) (See page 32.)

Question

I used the appliance after completing quick charging, but the charging indication/ notification lamp blinks immediately?

When the remaining capacity of the battery is low (charging indication/notification lamp blinks), does the flash power drop?

Can the appliance be used by men?

Dirt on the light guide cannot be removed.

Answer

If the appliance is used immediately after completing quick charging, the charging indication/notification lamp blinks. (Approximately 40 flashes possible) (See pages 32 and 33.)

Regardless of the remaining capacity of the battery, the flash power is the same.

Men can use the appliance apart from on the beard, but because generally men have more and thicker body hair than women, the stimulation of the skin may be stronger. If the dirt cannot be removed even after cleaning, contact

purchased the unit or a service center authorized by Panasonic for repair.
[Repair part]
Body attachment:

the store where you

WESWH93S3107

Troubleshooting **Problem** Action Check that the level indication lamp (blue) has changed from blinking to glowing. Press the flash switch perpendicularly against the The flash window skin until it flashes. (See does not flash. page 35.) If 10 continuous flashes have been emitted, temporarily move the appliance away from the skin The flash preparation time The flash preparation may become longer to time is long. prevent the main body from overheating. The power turned off Use the appliance in an ambient temperature after all lamps blinked and the appliance approximately 5 °C or higher. beeped twice. (See page 37.) For safety reasons, the power will turn off automatically if The power turned off the appliance is not used for approximately 1 minute after during use. the power is turned on. Turn the power on again.

Problem		Action
The appliance can only be used a few times after charging.	•	The battery has reached the end of its life. (See this page.)
The appliance does not operate even after pressing the power switch.	•	Press and hold the power switch (for approximately 3 seconds) when turning the power on.
	•	Charge the appliance. (See page 32.)

If the problems still cannot be solved, contact the store where you purchased the unit or a service center authorized by Panasonic for repair.

Battery life

Approximately 3 years.

If the operating time is significantly shorter even after a full charge, the battery has reached the end of its life. (The battery life may change due to usage and storage conditions.)

The battery in this appliance is not intended to be replaced by consumers. Have the battery replaced by an authorized service center.


Removing the built-in rechargeable batteries

Remove the built-in rechargeable batteries before disposing of the appliance.

Please make sure that the batteries are disposed of at an officially designated location if there is one.

This figure must only be used when disposing of the appliance, and must not be used to repair it. If you dismantle the appliance yourself, it may cause malfunction.

- Remove the appliance from the AC adaptor.
- If the appliance can be turned on and emit light, use it until the battery is depleted.
- Wait for at least 10 minutes after turning the power off before disassembling the appliance.
- Perform steps ① to ⑧ and lift the battery, and then remove it.
- Take care not to short circuit the positive and negative terminals of the removed battery, and insulate the terminals by applying tape to them.

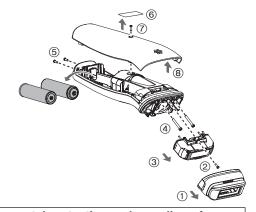For environmental protection and recycling of materials

This appliance contains Li-ion batteries.

Please make sure that the batteries are disposed of at an officially assigned location, if there is one in your country.

Specifications

Product name	IPL Hair Removal System
Product model	ES-WH93
Rated input	100-240 V∼ 50-60 Hz 180 mA
Rated output	4.8 V=== 1.25 A
Charging time	Approx. 3 hours
Place of origin	Main body: Japan AC adaptor: China
Date of production	See engraving on the commodity.
Execution standard	GB4706.1, GB4706.9, GB4706.15

This product is intended for household use only.

MEMO	

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