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Thank you for choosing a Panasonic Electric Toothbrush. Please read all instructions before use.

Safety Precautions

Warning

- Keep out of the reach of toddlers and infants. Failure to do so may result in accidents or accidental ingestion of the accessories or removable small parts.

Caution

- People who have severe periodontal disease should consult a doctor before use. Failure to do so may cause damage to teeth and gums. Do not use the brush in combination with toothpaste that contains a lot of abrasive, such as toothpaste for removing tar or for whitening. Do not force the bristles hard against the teeth or gums, or apply the bristles to one area for prolonged periods. Do not apply the solid part of the brush head to teeth. Toddlers, infants, people who are unable to operate and people with weak oral senses should not use this appliance. If you have oral health problems such as teeth undergoing treatment or unstable dentures, avoid using the brush on those parts as doing so may cause damage to teeth and gums. Never modify, disassemble, or repair the appliance yourself as doing so may cause unexpected accidents or trouble. Clean using a mild detergent. Never use thinners, benzene, or alcohol to clean the toothbrush since it may cause malfunction or cracking/discholoration of the parts. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.

Parts identification

- Main unit: Power switch [0/1], Battery cover, Dual-Edged Brush (x1), Inverted V-shaped bristles for interdental, Extra fine bristles for between teeth and gums. Cap.

Attaching and removing the cap. Attach and remove the cap straight.

Attaching the brush head

- Align the mark on the brush head with the mark on the main unit, and rotate it in the direction of the arrow until it clicks. Inserting the brush head without aligning the mark on the brush head with the mark on the main unit may damage the brush head or the main unit. Tightening it too strongly may cause deformation and damage of the main unit and the brush head.

Removing the brush head

- Rotate the brush head and align the marks. Pull to remove the brush. Pulling the brush head while it is fastened may damage the brush head or the main unit. When toothpaste gets hard and makes it difficult to remove the brush head, put it in water or warm water for a while and then remove.

Replacing the battery

When the battery is low, there is less vibration and a weaker brushed feeling. Please replace the alkaline battery with a new one.

- Remove the battery cover by rotating it. Insert an alkaline LR03 (AAA) alkaline battery. It is possible to use the appliance for approximately 90 days using a Panasonic alkaline LR03 (AAA) alkaline battery (using the appliance once a day for 2 minutes).

Attach the battery cover by rotating it.

- Batteries are not supplied. Please purchase separately. Before replacing the alkaline battery, ensure the unit is thoroughly dry to prevent water drops entering the main unit. Always use alkaline batteries. Use of non-alkaline batteries may result in a drastically reduced operating time. Do not use rechargeable batteries. We recommend that you use a Panasonic LR03 (AAA) alkaline battery. Always follow the cautions printed on the battery. Check the polarities (+ and -) of the battery and install it correctly. Remove the battery when not using the appliance for an extended period of time. Please dispose of used battery appropriately. Do not disassemble the battery or throw it into a fire. Do not short-circuit the battery.

Using the toothbrush

- Remove the cap. Confirm that the brush is fixed firmly.

Replacing the battery

- Place the toothbrush in your mouth. Press the power switch firmly to operate the toothbrush. Press the power switch with the ball of the finger. Do not press with a fingernail. Doing so may tear the rubber of the power switch. Brush your teeth by moving the main unit by hand. Brush by moving the toothbrush in the same way as when you brush by hand. Just applying the toothbrush to teeth cannot fully remove plaque, so it is recommended that you move the toothbrush. Do not apply the brush head to teeth or gums too strongly. When finished brushing, stop the motion by pressing the power switch firmly. In order to prevent splattering of saliva, etc., stop operation of the toothbrush while it is still in the mouth.

Tips on usage

Boundary between the tooth surfaces or teeth and gums

- Press the bristles at a 45° angle. Place the brush so that the minute side bristles can fit in between the teeth and gums.

At the back of the front teeth or misaligned teeth

- Press the brush vertically. Place the brush according to the angle of the teeth.

After use

Washing the brush

- Wash the brush under running water. Do not wash with water hotter than 80 °C. Do not pull the bristles of the brush. (Doing so may cause bristles to come out.) Shake well to remove water after washing.

Washing the main unit

- Make sure that the battery cover is closed. Wash the main unit under running water. Do not wash with water hotter than 40 °C. Do not immerse the main unit under water. (Otherwise water may get inside.)

Washing the cap

- Wash the Cap under running water. Do not wash with water hotter than 40 °C.

Changing the brush head (Sold separately)

Replacing the brush every 3 months is recommended. If the bristles are bent to one side or worn, the toothbrush is not very effective for brushing.

Table with 3 columns: Part number, Replacement brush head for EW-DS11, Quantity. Row 1: WEW0957-W (White), [Image of brush head], 2 pcs.

Specifications

Table with 2 columns: Product name, Electric toothbrush. Product model, EW-DS11. Rated voltage, 1.5V---. Battery, LR03 (AAA) alkaline battery x 1 (sold separately). Weight, Approx. 35 g (including brush and cap) (excluding battery). Place of origin, China. Date of production, See engraving on the commodity. Execution standard, Company standard. Q/SXDGWB 4.

Notes

- Do not use for any purpose other than dental care. Do not drop or handle roughly. Please be careful not to get hair products, hand cream, etc. on the main unit because it may ruin the appearance.

Tips on usage

Bite (back teeth)

- Apply the brush to the teeth at an angle of 90° from the top.

Notes

- Do not use for any purpose other than dental care. Do not drop or handle roughly. Please be careful not to get hair products, hand cream, etc. on the main unit because it may ruin the appearance.

After use

Washing the brush

- Wash the brush under running water. Do not wash with water hotter than 80 °C. Do not pull the bristles of the brush. (Doing so may cause bristles to come out.) Shake well to remove water after washing.

Washing the main unit

- Make sure that the battery cover is closed. Wash the main unit under running water. Do not wash with water hotter than 40 °C. Do not immerse the main unit under water. (Otherwise water may get inside.)

Washing the cap

- Wash the Cap under running water. Do not wash with water hotter than 40 °C.

Tips on usage

Bite (back teeth)

- Apply the brush to the teeth at an angle of 90° from the top.

Notes

- Do not use for any purpose other than dental care. Do not drop or handle roughly. Please be careful not to get hair products, hand cream, etc. on the main unit because it may ruin the appearance.

After use

Washing the brush

- Wash the brush under running water. Do not wash with water hotter than 80 °C. Do not pull the bristles of the brush. (Doing so may cause bristles to come out.) Shake well to remove water after washing.

Washing the main unit

- Make sure that the battery cover is closed. Wash the main unit under running water. Do not wash with water hotter than 40 °C. Do not immerse the main unit under water. (Otherwise water may get inside.)

Washing the cap

- Wash the Cap under running water. Do not wash with water hotter than 40 °C.

Note

- Shake the cap and brush head well to drain off water, and dry using a dry towel, etc. After this, put the cap on the main unit. Leaving them wet may dirty or damage the pouch or bag.

Changing the brush head (Sold separately)

Replacing the brush every 3 months is recommended. If the bristles are bent to one side or worn, the toothbrush is not very effective for brushing.

Table with 3 columns: Part number, Replacement brush head for EW-DS11, Quantity. Row 1: WEW0957-W (White), [Image of brush head], 2 pcs.

Specifications

Table with 2 columns: Product name, Electric toothbrush. Product model, EW-DS11. Rated voltage, 1.5V---. Battery, LR03 (AAA) alkaline battery x 1 (sold separately). Weight, Approx. 35 g (including brush and cap) (excluding battery). Place of origin, China. Date of production, See engraving on the commodity. Execution standard, Company standard. Q/SXDGWB 4.

Notes

- Do not use for any purpose other than dental care. Do not drop or handle roughly. Please be careful not to get hair products, hand cream, etc. on the main unit because it may ruin the appearance.

After use

Washing the brush

- Wash the brush under running water. Do not wash with water hotter than 80 °C. Do not pull the bristles of the brush. (Doing so may cause bristles to come out.) Shake well to remove water after washing.

Washing the main unit

- Make sure that the battery cover is closed. Wash the main unit under running water. Do not wash with water hotter than 40 °C. Do not immerse the main unit under water. (Otherwise water may get inside.)

Washing the cap

- Wash the Cap under running water. Do not wash with water hotter than 40 °C.

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English

Using the toothbrush

- Remove the cap. Confirm that the brush is fixed firmly.

Replacing the battery

- Place the toothbrush in your mouth. Press the power switch firmly to operate the toothbrush. Press the power switch with the ball of the finger. Do not press with a fingernail. Doing so may tear the rubber of the power switch. Brush your teeth by moving the main unit by hand. Brush by moving the toothbrush in the same way as when you brush by hand. Just applying the toothbrush to teeth cannot fully remove plaque, so it is recommended that you move the toothbrush. Do not apply the brush head to teeth or gums too strongly. When finished brushing, stop the motion by pressing the power switch firmly. In order to prevent splattering of saliva, etc., stop operation of the toothbrush while it is still in the mouth.

Tips on usage

Boundary between the tooth surfaces or teeth and gums

- Press the bristles at a 45° angle. Place the brush so that the minute side bristles can fit in between the teeth and gums.

At the back of the front teeth or misaligned teeth

- Press the brush vertically. Place the brush according to the angle of the teeth.

After use

Washing the brush

- Wash the brush under running water. Do not wash with water hotter than 80 °C. Do not pull the bristles of the brush. (Doing so may cause bristles to come out.) Shake well to remove water after washing.

Washing the main unit

- Make sure that the battery cover is closed. Wash the main unit under running water. Do not wash with water hotter than 40 °C. Do not immerse the main unit under water. (Otherwise water may get inside.)

Washing the cap

- Wash the Cap under running water. Do not wash with water hotter than 40 °C.

Note

- Shake the cap and brush head well to drain off water, and dry using a dry towel, etc. After this, put the cap on the main unit. Leaving them wet may dirty or damage the pouch or bag.

Changing the brush head (Sold separately)

Replacing the brush every 3 months is recommended. If the bristles are bent to one side or worn, the toothbrush is not very effective for brushing.

Table with 3 columns: Part number, Replacement brush head for EW-DS11, Quantity. Row 1: WEW0957-W (White), [Image of brush head], 2 pcs.

Specifications

Table with 2 columns: Product name, Electric toothbrush. Product model, EW-DS11. Rated voltage, 1.5V---. Battery, LR03 (AAA) alkaline battery x 1 (sold separately). Weight, Approx. 35 g (including brush and cap) (excluding battery). Place of origin, China. Date of production, See engraving on the commodity. Execution standard, Company standard. Q/SXDGWB 4.

Notes


- Do not use for any purpose other than dental care. Do not drop or handle roughly. Please be careful not to get hair products, hand cream, etc. on the main unit because it may ruin the appearance.

使用说明书 Operating Instructions

(家用) 电动牙刷 (Household) Electric Toothbrush

型号 EW-DS11 Model No.

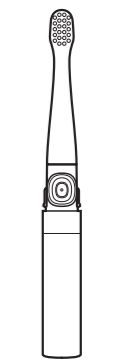中文 ..... 2 English ..... 8

在使用本产品前, 请仔细阅读本说明书并妥善保管, 以备将来参考。 Before operating this unit, please read these instructions completely and save them for future use.

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© CH, EN