# Panasonic

## Operating Instructions

Sonic Vibration Toothbrush

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

### EW1035

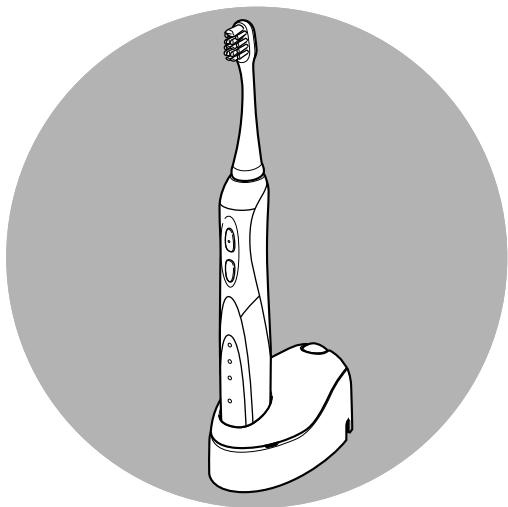
Before operating this product, please read these instructions completely and save this manual for future use.

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










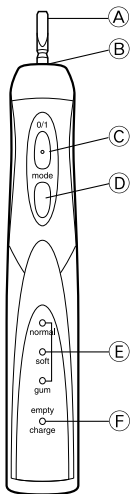
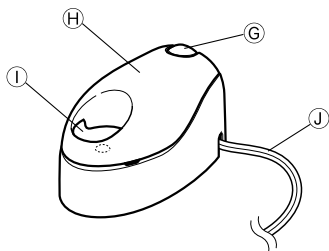
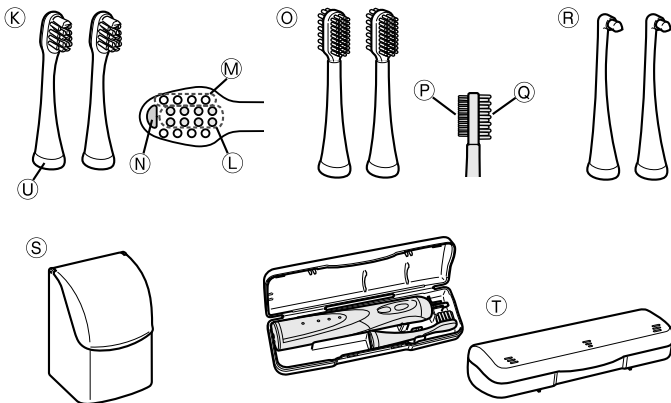
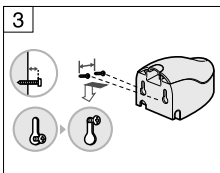
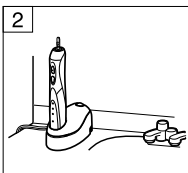
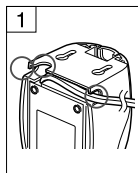
ENGLISH — Page 4

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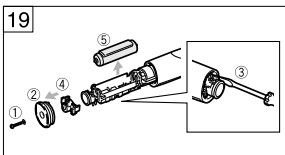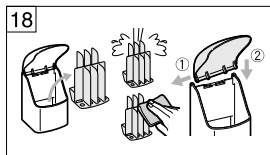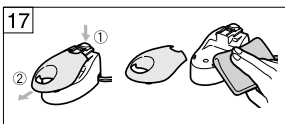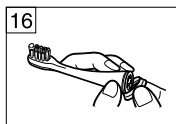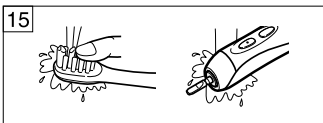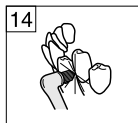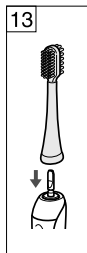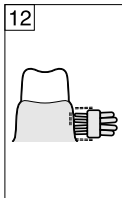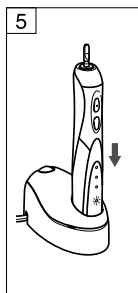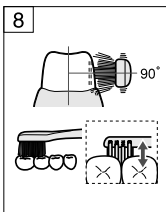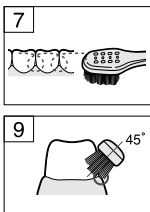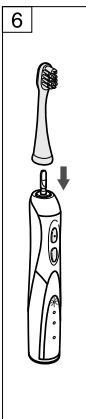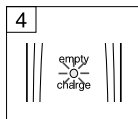## Safety Instructions

Please read the safety instructions before using this product. The safety instructions provided here are to help you use this product safely and correctly, and also to prevent you from injuring yourself or others.

The safety instructions shown here are classified as follows.

**WARNING:** Indicates warnings that, if not heeded, may lead to death or serious injury.

**CAUTION:** Indicates precautions that, if not followed, may lead to serious injury or property damage.

### Illustrated indications

⊘ indicates actions that you must NOT take.

● indicates actions that you should take or instructions that you should follow.

### WARNING

Unplug

- Be sure to switch off the product and remove the plug from the electrical outlet before performing any maintenance. Do not touch or pull out the plug with wet hands. Failure to do so may result in electrical shock or other injury.

Be sure to follow the instructions below.

- Clean the plug periodically so that dust does not accumulate in it. Such dust buildup may trap

moisture over time and reduce the effectiveness of the insulation, thus resulting in overheating and fire.

- Use the product at the rated voltage. Failure to do so may result in accidents such as fire, electrocution and/or personal injury, including burns.

Prohibited use

- Never operate this product if it has a damaged cord or plug, if the cord or plug is warm, or if the plug fits loosely in the electrical outlet. Failure to do so may result in electrocution, short-circuit and/or fire.
- Do not damage, break or modify the cord. Do not bend, pull or twist the cord with unnecessary force. Do not place heavy objects on the cord or allow the cord to get pinched between objects. Doing so may cause fire, electrocution, short-circuit and/or burns.

Never disassemble the product.

- Do not modify, disassemble or repair the product at any time except when disposing of it. Doing so may result in fire, electrocution and/or personal injury. For repair or maintenance, take the product to the store where you purchased it or to an authorized service centre.

Do not use in areas where water is present.

- Do not use the charger in the bath. Doing so may result in electrocution and/or short-circuit.
- Do not wash the charger. Doing so may result in electrocution.
- Do not reach for the charger that has fallen into water. Unplug the charger immediately to avoid the risk of electrocution.

## CAUTION

Be sure to follow the instructions below.

- Use this product only for its intended use as described in this manual. Failure to do so may result in accidents and/or personal injury.
- Keep the product out of reach of infants. Failure to do so may result in burns and/or electrocution.
- Use the accessories recommended by the manufacturer. Failure to do so may result in accidents and/or personal injury.

Unplug

- Unplug when the product will not be used for more than a week. Failure to do so may result in electrocution and/or fire caused by electrical leakage as a result of inadequate insulation.

Prohibited use

- Never unplug it by pulling the cord. Doing so may result in electrocution and/or short-circuit.
- Discontinue use immediately if the product has been dropped, hit or otherwise damaged through impact. Such damage may cause electrocution and/or fire.
- This product should not be used by infants, or by those who have limited sensation inside their mouth. Doing so may result in accidents or personal injury.
- Don't let infants use. Doing so may result in personal injury.
- Do not place the toothbrush forcibly against your teeth or gums, or do not leave the toothbrush in one location in the mouth for a long time. Doing so may damage the teeth and gums. People with tooth or gum disease or people receiving oral treatment should consult their dentist before use.
- Do not brush loose dentures or teeth that are under dental treatment. Also, keep the brush away from any sores or other sensitive areas inside the mouth. Failure to observe these precautions may aggravate the condition and/or result in an accident.
- Do not use this toothbrush with toothpaste that contains a large amount of abrasives for removing stains from teeth, or with toothpaste that contains a lot of whitening agents. Failure to observe these precautions may damage your teeth or gums.

- Do not place any metal objects, such as coins and clips, on the charger. Doing so may cause the objects to heat up, resulting in burns.
- Do not place the plastic portion of the toothbrush against the teeth or gums. Doing so may damage your teeth or gums.

## For your personal dental care

Thank you for purchasing the Sonic Vibration Toothbrush.

The Sonic Vibration Toothbrush is a new kind of electric toothbrush designed for those who wish to keep their teeth clean and healthy. This product will enable you to maintain good dental hygiene, which is as important as other forms of personal hygiene.

### Sonic Vibration Toothbrush Technology

The Sonic Vibration Toothbrush incorporates innovative “linear motor” technology that uses magnetic attraction and repulsion to provide outstanding performance. The brush’s 31,000 horizontal and 34000 vertical strokes per minute help remove plaque gently yet effectively.

### Other benefits

The toothbrush stops brushing briefly after 2 minutes of use and starts again as a convenient reminder of the recommended brushing time.

The intensity of the toothbrush vibration is automatically weakened when the brush is pressed too hard against your teeth.

## PARTS IDENTIFICATION

- Ⓐ **Toothbrush shaft**
- Ⓑ **Waterproof rubber seal**
- Ⓒ **On/off button**
- Ⓓ **Mode selector**

- Ⓔ **Mode indicator**
- Ⓕ **Charging indicator**
- Ⓖ **Charger cover release button**
- Ⓗ **Charger cover**
- Ⓘ **Charging section**
- Ⓙ **Power cord**
- Ⓚ **Multi-reach brush**
- Ⓛ **Float mechanism**
- Ⓜ **Minute side bristles**
- Ⓝ **Dense tip bristles**
- Ⓞ **Double-sided silicone brush**
- Ⓟ **Normal**
- Ⓠ **Hard**
- Ⓡ **Precision-tip brush**
- Ⓢ **Brush storage case**
- Ⓣ **Carrying case**
- Ⓤ **Identification ring**

## BEFORE USE

This high-speed, vibrating, electric toothbrush is different from conventional electric toothbrushes. Therefore:

- You only need to lightly press the bristles against your teeth for brushing.
- When first using this product, some may experience a tickling or tingling sensation caused by the high-speed vibrations. The sensation should subside after you use the toothbrush several more times.

## Installing the charger

### Cord (Fig. 1)

Place the cord in the right, left or center cord guide, depending on where the charger will be placed.

### **Installing near a sink or on a shelf (Fig. 2)**

Place the charger on a level, secure surface. Insert the plug into an electrical outlet.

### **Mounting on the wall (Fig. 3)**

1. Purchase two round wood screws (diameter: 2.7 – 3.1 mm).

2. Decide where you want to mount the charger.

Fasten the screws so that they extend about 2.5 mm from the wall and are positioned 28 mm apart.

3. Align the holes on the back side of the charger with the screws and secure in place.

\*Mount the brush storage case in the same manner, if desired.

### **CAUTION:**

- Do not use the charger in a shower stall or other area with running water that come in direct contact with the charger.
- Do not remove the plug from the electrical outlet with wet hands.
- Be sure to use only a 100-240 VAC power supply. (110-120 VAC for Taiwan)

### **NOTE:**

- When mounting the charger on the wall, be sure to insert the screws in the wall studs or other sturdy section of the wall.

## **CHARGING**

Be sure to charge the main unit for at least 15 hours:

- 1) before using for the first time;
- 2) when you have not used the toothbrush for more than three months; or

3) if the charging indicator blinks green during use (Fig. 4)

\*The toothbrush can be used one or two times even after the charging indicator blinks green depending on usage conditions.

\*When the charging indicator blinks green, the mode indicator does not light up. (Mode selection is still possible even though the mode indicator does not light.)

Be sure to turn off the toothbrush before charging.

Place the toothbrush upright in the charger. (Fig. 5)

The charging indicator will light up red.

- The toothbrush charges whenever it is placed in the charger.
- The indicator light does not go off even after the toothbrush is fully charged.
- The toothbrush will provide about 40 minutes of use when the battery is fully charged.
- The operating duration will become shorter as the battery gets older.

### **CAUTION:**

Do not leave metal objects, such as coins or paper clips, on top of the charger. They may become very warm.

### **NOTE:**

- During charging, the toothbrush will get warm, but this is normal.
- Use only the exclusive charger provided.
- If the toothbrush operates only for a short time after being fully charged, this means that the battery is nearing the end of its life cycle. Take the unit to an authorised service centre to have the

battery replaced.

- Use the toothbrush at an ambient room temperature of 0° to 40°C.

## HOW TO USE

1. Place the brush head firmly onto the toothbrush shaft. (Fig. 6)
  - There is a small space between the brush and the main unit. This is normal and is necessary for the head to vibrate properly.
2. Place the toothbrush in your mouth.
  - Apply small amount of normal toothpaste to the toothbrush, if desired, before turning the power on to prevent the toothpaste from spattering.
3. Press the on/off button to turn on the toothbrush.
  - When the toothbrush first starts, the “soft-start” function causes it to take about 2 seconds to reach its normal speed.
4. Press the mode selector to select the mode.
  - Pressing the mode selector toggles the operation mode in the order of “normal” → “soft” → “gum”.

### Mode: normal

Use: Standard brushing

Brush movement:

Combined movement of horizontal, vertical and diagonal.

### Mode: soft

Use:

- When teeth or gums are sensitive.
- When brushing in the Normal mode is too strong.

Brush movement:

Horizontal movement only.

### Mode: gum

Use: For gums using the double-sided silicone brush.

Brush movement:

Horizontal movement only.

Note:

You can use the double-sided silicone brush in the normal and soft modes, too.

5. When you have finished brushing, press the on/off button to turn off the power.
  - To prevent foam from spattering, turn the toothbrush off before removing it from your mouth.

### Brushing duration notification timer

- The toothbrush stops temporarily after 2 minutes of use, and starts again. This reminds you that you should brush your teeth for about 2 minutes.
- When the mode selector is pressed, the toothbrush provides 2 minutes of timed operation. In the gum mode, the timer does not operate.

### Strong pressure preventive function

Pressing the brush bristles too hard against your teeth will activate this feature, causing the intensity of the vibration to be weakened automatically as a safeguard against excessive pressure.

### PRECAUTION:

- When you first begin using the toothbrush, you may experience slight bleeding of the gums even though your gums may be

healthy. This is because your gums are receiving stimulation by the toothbrush for the first time. The bleeding should stop after one or two weeks of use. Bleeding that persists for more than two weeks may be an indication of an abnormal condition of your gums, such as pyorrhea. In this case, stop using the toothbrush and consult your dentist.

- Do not allow hair products (e.g., gel or oil) or hand lotion to come in contact with the unit as they may damage its surface.

#### **NOTE:**

Replacing the brush head every 3 months is recommended.

If the bristles are bent to one side, the toothbrush is not very effective for brushing. Replace the brush with a new one when the bristles are bent or worn.

## **TIPS ON USAGE**

### **Multi-reach brush (EW0923)**

For brushing your teeth.

Place the bristles of the toothbrush gently against your teeth and move the toothbrush slowly in a circular motion (Fig. 7).

1. For the side teeth and back teeth (Fig. 8)
  - Place the bristles at a 90° angle.
  - Place the brush so that the bristles of the float mechanism, which features miniature springs in the base to provide fluid brushing action, can fit in

between the teeth.

2. For between teeth and gums (Fig. 9)
  - Place the bristles at a 45° angle.
  - Place the brush so that the minute side bristles can fit in between the teeth and gums.
3. For reaching the far side of the back teeth (Fig. 10)
  - Place the brush so that the dense tip bristles can reach the molars.
4. For the back side of front teeth or where teeth are not straight (Fig. 11)
  - Brush vertically.
  - Move toothbrush to match the angle of the teeth.

### **Double-Sided Silicone Brush (EW0931) (Fig. 13)**

The double-sided silicone brush is for your gums.

- Place the brush lightly against your gums.
- Place the bristles at a 90° angle. (Fig. 12)
- The gum mode is recommended for use with the silicone brush. In the gum mode, the brush movement changes automatically. (The silicone brush can be used in the normal and soft modes, too.)
- Do not place the silicone brush onto your teeth. Doing so may damage the bristles of the silicone brush.

### **Using the double-sided silicone brush**

- Attach the brush so that the normal side faces the front side of the unit. The brush cannot be

attached with the hard side facing the front.

- Use the desired side of the brush: normal or hard.
- The hardness of the bristles of the normal side of the silicone brush is the same as that of the optionally available silicone brush EW0930.

### **Placing the brush**

- Be sure to place the tips of the bristles so they are lightly touching your gums. (Pressing the bristle tips hard against your gums may stop the motion of the brush.)
- You will not feel the brushing action as much as you do when using a nylon-bristle brush.

### **Brushing duration**

- Approximately two minutes of use once a day is recommended.
- Use the silicone brush for 10 to 20 seconds on sections of the gums that are of special concern.
- For other gum areas, move the brush slowly and place the entire gum area.

### **Gum mode is recommended for those who:**

Are experiencing bleeding of the gums;  
Are experiencing gum discomfort or pain periodically;  
Wish to improve their oral hygiene and freshen their breath.

### **NOTE:**

Do not pull or bite hard on the silicone brush. Doing so may damage the bristles.

### **Precision-tip brush (EW0940) (Fig. 14)**

For where teeth are not straight or for teeth that cannot be brushed adequately using the multi-reach brush.

Place the brush lightly against your teeth.

- Place the brush vertically.
- Move slowly along between the teeth.

### **AFTER USE (Fig. 15)**

#### **Toothbrush**

Wash the brush under running water by moving the float (green bristles) section.

- Shake well to remove water after washing.
- Do not pull the bristles of the multi-reach brush. (Pulling the bristles may dislodge them, causing them to fall out.)
- Do not spread apart the bristles of the double-sided silicone brush in order to wash it. (Doing so may damage the silicone brush.)
- Do not shake the main unit strongly with the brush attached. (Doing so may damage the main unit and the brush.)

#### **Main unit**

Wash the main unit under running water.

- Do not immerse the main unit in water.
- Do not rub or damage the waterproof rubber seal. (Doing so may cause water to enter the main unit.)
- Do not wash with water hotter than 40°C.

After washing with water, wipe dry with a cloth. Return the main unit to the charger. Return the brush to the brush storage box.

If the identification ring comes off, put it back on. (Fig. 16)

- Do not forcibly pull the identification ring since it can break easily.
- Using the toothbrush without the identification ring has no effect on performance; the ring is simply for user identification purposes.

- Do not use harsh chemicals such as thinner, benzene or alcohol as these will discolour the surface.

#### **CAUTION:**

- Do not wash or submerge the charger in water.
- Clean the plug periodically. Dirt and dust buildup may cause electrical shock or reduce the effectiveness of the insulation, thus resulting in overheating and a fire.

## **MAINTENANCE**

#### **CAUTION:**

Be sure to remove the plug from the electrical outlet before cleaning the product. Failure to do so may result in electrical shock or other injury.

#### **Charger (Fig. 17)**

1. Pressing the charger cover release button, pull the charger cover to remove.
2. Using a soft cloth, wipe off any dirt on the charger and the charger cover.
3. About once every six months wipe the prongs of the plug with a dry cloth.

#### **Brush storage box (Fig. 18)**

1. Remove the separator.
  2. Wash the separator with water and dry with a cloth.
- If the lid comes off, ① attach one of the hooks first and then ② attach the other hook from the back of the storage case.
  - If cleaning requires a detergent, use a mild, neutral cleanser.

## TROUBLESHOOTING

<b>Problem</b>	<b>Possible cause</b>	<b>Remedy</b>
The intensity of the vibration weakens.	The toothbrush is pressed hard against your teeth and excessive pressure preventive function is activated.	Place the brush gently against your teeth.
After turning on the power, the toothbrush vibration speed increases gradually.	This is normal operation of the soft-start function. (See page 9.)	Not a malfunction. Continue using.
After 2 minutes of use, the toothbrush stops temporarily and starts again.	The timer is reminding you that you have been brushing for 2 minutes. (See page 9.)	Not a malfunction. Continue using.
The toothbrush operates only for a short period of time even immediately after being charged.	The main unit was not positioned straight up in the charger and was unable to charge properly. (See page 8.)	Stand the main unit up straight so that the bottom touches the charging base, and make sure the charging indicator comes on.
	The battery has reached the end of its life cycle (about 3 years).	Take the unit to an authorised service centre to have the battery replaced. (Do not attempt to replace the battery yourself since you may damage the unit's waterproofing or internal mechanisms.)
	Charging time was insufficient.	Charge the main unit for at least 15 hours. (See page 8.)

<b>Problem</b>	<b>Possible cause</b>	<b>Remedy</b>
The toothbrush doesn't operate.	You have just purchased the product or it has not been used for more than 3 months.	Charge the main unit for at least 15 hours. (See page 8.)
The direction of the brush movement changes. (The brush head moves from side to side.)		Not a malfunction. Remove the brush head once and put it back to the correct position.
The toothbrush is not working properly. (For example, its movement is uneven.)		Place the main unit in the charger for at least 15 hours, then turn it on again to see if it operates properly.

If the product still does not function properly, stop using it and contact an authorised service centre.

**WARNING:**

Do not disassemble the product at any time other than when disposing of it. Doing so may result in short circuit, fire, or injury.

For repair or maintenance (such as battery replacement), take the product to an authorised service centre.

The supply cord of the charging adaptor cannot be replaced; if the cord is damaged, the entire charger should be discarded.

## DISPOSAL OF THE PRODUCT (Fig. 19)

The power source of this product is nickel-metal hydride battery which can be recycled. Please dispose of the device and the battery in accordance with national or public regulations or environmental protection or recycling regulations.

### WARNING:

Do not disassemble the product at any time other than when disposing of it. Failure to observe this point may result in short circuit, fire, or injury.

For repair or maintenance (such as battery replacement), take the product to an authorised service centre.

Removing the Ni-MH battery before disposing of the product

- Before removing the battery, be sure to completely drain the unit of power.
- Use a screwdriver to disassemble the unit to dispose of the battery. Remove ① and ②. Pressing the toothbrush shaft (see page 2), remove the two hooks ③. Remove the ④ and ⑤ as shown in Fig. 19.
- Be careful not to short-circuit the negative and positive poles of the battery.

## SPECIFICATIONS

Rated voltage:

100 - 240 V AC  
(110-120 VAC for Taiwan)  
50 - 60 Hz

Battery:

1 Ni-MH battery (1.2 V DC; 1,000 mAh)

Charging time: 15 hours

Input current: 35 mA

Power consumption:

2 W (when fully charged)

Brush rotating speed:

approx. 31,000 horizontal strokes/min.  
approx. 34,000 vertical strokes/min.

Brush movement (with no load):

normal mode: approx. 1.0 mm horizontal (width) and approx. 6° vertical (height) movement  
soft mode: approx. 0.7 mm horizontal (width) movement  
gum mode: approx. 0.7 – 1.00 mm (changes automatically)

Operating time:

approx. 40 minutes of continuous use (with 300 g load at a room temperature of 20°C after the battery has been fully charged.)

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**Panasonic Electric Works Co., Ltd.**